

iMassage

Smart Eye Massager


iGear
innovation. delight



**Rejuvenate your eyes
anytime anywhere !**

Model No.: iG 1047

EFFECTIVE FOR:



DRY EYES
TIRED EYES



STRESS
ANXIETY



MIGRAINES
HEADACHES

THE BENEFITS

The iMassage is effective against:

Migraines, Sinus Pressure & Headaches

Dry & tired eyes

Sleeping problems

Stress & Anxiety

Eye bags & dark under-eye circles

BEFORE



AFTER



RESULTS

Experience instant and long-term benefits:

No Migraines or Headaches

Healthy & fresh eyes

Improves sleep quality

No Stress or Anxiety

No eye bags & dark under-eye circles



THE SCIENCE BEHIND iGear iMASSAGE

HOW DOES IT WORK?

The iGear iMassage Smart Eye Massager uses safe air compression technology and adopts kneading, trigger point therapy and rhythmic percussion massaging to give you the gentle touch you need to battle against headaches, eye strain, sleeping problems, stress and more.

The built-in heating pads provide a soothing adjustable temperature between 100 to 107 degrees Fahrenheit (38°C - 42°C) which helps to stimulate blood circulation around your eyes, reducing eye puffiness, dark circles, improving your sleep quality and making your eyes feel healthy and refreshed!



WHO IS IT FOR?

If you suffer from headaches, dry and tired eyes, bad quality sleep or stress then the iGear iMessage Smart Eye Massager is perfect for you.

It's been designed to keep your eyes healthy and relieve all eye-related problems including migraines, sinus pressure, eye fatigue, puffiness, eye bags and more, to be used in the comforts of your own home (or wherever you prefer) and avoid expensive migraine medicine purchases & costly massage therapist visits.

It's the perfect solution for anyone wanting a life with healthy eyes, great sleep & no headaches!

FAQ

1. How long does it take to feel the relief?

The iGear iMessage Smart Eye Massager provides both immediate relief (in as little as 15 minutes a day) as well as long-term results with repeated use.

2. How well does this work on migraines?

Wonderfully. The heating function of the massager is a true saviour against migraines.

3. How frequent should I be using this device per week?

We recommend to use the massager daily or every other day for lasting relief.

4. Is the iGear iMessage Smart Eye Massager Safe?

The device is 100% safe to use! Our device has received all necessary certifications to be considered as a completely safe device for you and your eyes.

USAGE METHOD

Step 1

Press the power button for 3 seconds to start the Eye massager. Once On, the massager will by default start in the 1st massager mode i.e., Heat compress + Vibration + Airbag. Press the power button once to switch from the massager mode (Level 1 heat compress + vibration + airbag) to heat mode (Level 1 heat compress only)

Step 2



Press the power button once to switch to 2nd massager mode (Level 2 heat compress + vibration + airbag). Press the power button once to switch to Level 2 heat mode (Level 2 heat compress only)

Step 3

Press the power button once to switch to 3rd massager mode (Level 3 heat compress + vibration + airbag). Press the power button once to switch to Level 3 heat mode (Level 3 heat compress only)

Note: To switch off the massager, press the power button for 3 seconds and the massager will turn off.

iGear iMessage Smart Eye Massager Diagram

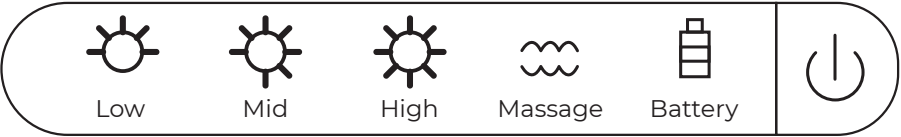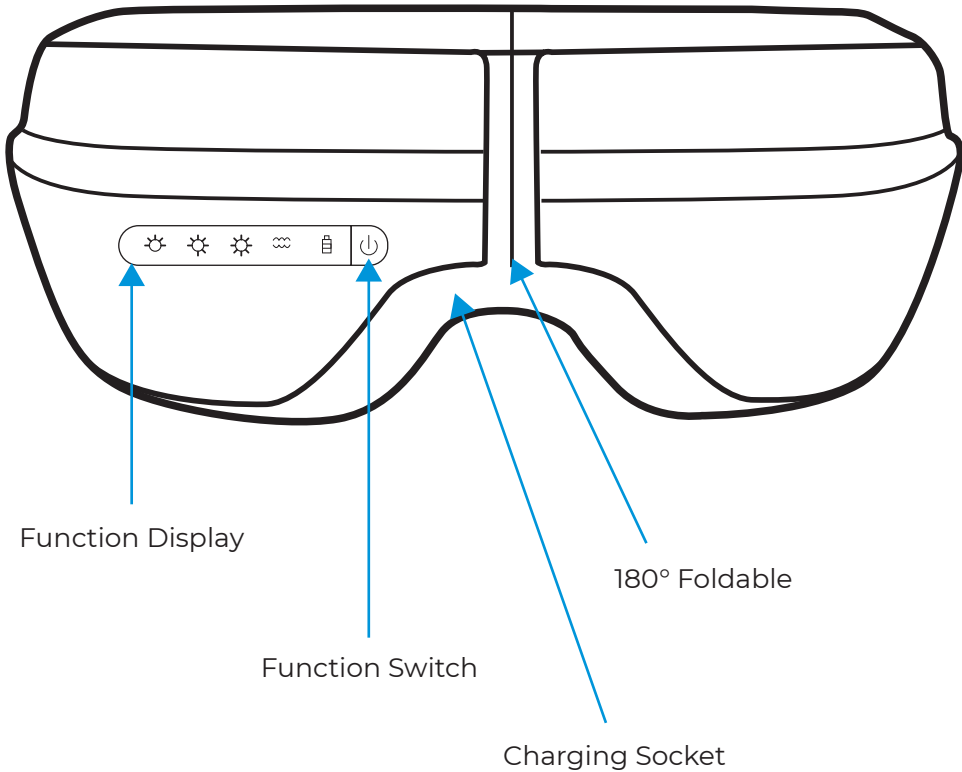



Matters Needing Attention

Do not use it for infants or people without the ability to express consciousness.

Do not use if there is eye injury or operation.

Do not use it when you have symptoms of glaucoma, cataract, retinal detachment and other eye diseases.

Do not use for stroke, hyper tension and hypotension.

Do not wear contact lenses when using this product.

It can not be used in high temperature, flammable, electro magnetic radiation and humid places.

Do not use when bathing or sweating.

The product has a heating surface, which should be paid attention to by people who are sensitive to heat.

Do not use with other medical electronic instruments, such as pace maker, artificial heart and lung life support medical electronic instruments.

The battery must be removed from the product before it is discarded. When removing the battery, the product must be powered off and the battery should be disposed of safely.

PACKAGE INCLUDES

- 1) iGear iMessage Smart Eye Massager
- 2) Instruction Manual
- 3) USB charging cable



Imported and Marketed By:
iGear Tech Pvt. Ltd.
204, Sai Corporate Park,
Laxmi Nagar, Off Link Road,
Goregaon West, Mumbai - 400104.

In case of consumer complaints:
Toll free number - 1800 22 8495
Email: sales@igear.asia
(Between 10am to 5pm, Mon- Fri)